INTELLIGENT
BIOLOGY™

TRI POLLAR

USER MANUAL

ENGLISH

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STOP MIRACLE

YOUNG

ENGLISH

**CONGRATULATIONS ON PURCHASING
YOUR TRIPOLLAR STOP MIRACLE YOUNG DEVICE!**

This user manual describes the TriPollar STOP MIRACLE YOUNG device and how to use it.

Supplementary information to be found at:
www.tripollar.com | cn.pollogen.com

TABLE OF CONTENTS

| | |
|--|-----------|
| ABOUT TRIPOLLAR RF TECHNOLOGY | 9 |
| TRIPOLLAR STOP MIRACLE YOUNG KIT INCLUDES | 9 |
| OVERVIEW OF THE TRIPOLLAR STOP MIRACLE YOUNG DEVICE | 10 |
| TREATMENT PROCEDURE | 11 |
| TROUBLESHOOTING | 13 |
| TECHNICAL SPECIFICATIONS | 14 |
| WARRANTY AND CUSTOMER SUPPORT | 15 |
| APPENDIX | 17 |

IMPORTANT NOTICE

© Pollogen Ltd. All rights reserved.

Design and technical specifications are subject to change without prior notice. Manufactured by Pollogen Ltd. Patented technology.

The names Pollogen, TriPollar are registered trademarks of Pollogen Ltd.

All intellectual property rights in this publication are owned by Pollogen Ltd. and protected by applicable copyright laws and international treaty provisions. Pollogen Ltd. retains all rights not expressly granted. No part of this publication may be reproduced in any form whatsoever or used to make any derivative work without prior written approval by Pollogen Ltd.

All other trademarks are the property of their respective owners. Other company and brand products and service names are trademarks or registered trademarks of their respective holders.

SAFETY

Important: Please read the full text of the user manual carefully prior to use of the product. Always adhere to all warnings and cautions. Be sure to follow all safety precautions below when using the TriPollar STOP MIRACLE YOUNG product.

To learn more about TriPollar RF technology, visit our website at:

www.tripollar.com | cn.pollogen.com.

- Only use this device for its intended purpose, as described in this manual and correctly use this device by strictly following the instructions of use. Failure to follow the instructions in this manual may cause damage or burns to your skin.
- Do not use the device over eyelids, closed eyes, thyroid area, lips, ears, breasts or genital areas.

- The sensation while using the device should be warm and pleasant. If you experience discomfort, discontinue use immediately.
- The use of this product may cause discomfort or adverse effects on people with allergic constitution, people with skin diseases and allergic dermatitis, people with sensitive skin, etc.
If adverse effects occur during use, you should stop using this appliance immediately.
- In case the device or the power supply generates excessive heat disconnect power supply from the electrical outlet & consult service.
- TriPollar STOP MIRACLE YOUNG is an electrical device that must be used properly.
To reduce the risk of injury from electric shock, adhere to all warnings and cautions below:

 **WARNINGS:**

- Do not touch the device should it fall into water.
- Do not place or store the device near any source of water, such as a sink or bathtub.
- Do not submerge the device under water.
- When plugged in, do not leave the device unattended.
- Before use, ensure that the system and its USB cord do not show any visible damage. Do not operate if the electric cord is damaged or if the system is disassembled.
- Only use the device with the electric cord provided.
- Keep the device out of the reach of children.
- Be sure that your hands are completely dry when operating the device.
- Turn the device OFF after use or before cleaning and maintenance.

CAUTIONS

- Use the device for personal use only.
- Use TriPollar Preparation Balm during treatment according to its instructions for use.
- Do not use TriPollar Preparation Balm after its expiration date.

DO NOT USE THE DEVICE IF ANY OF THE FOLLOWING APPLY TO YOU:

- Are under 18 years of age.
- Have a pacemaker or internal defibrillator.
- Have metal implants in the treatment area, such as gold threads. This does not apply to dental fillings or implants, which present no problem.
- Currently have cancer or a history of cancer, especially skin cancer or premalignant moles.
- Have severe concurrent conditions, such as cardiac disorders.
- Are pregnant or nursing.
- Have an impaired immune system due to immunosuppressive diseases such as AIDS, or are using immunosuppressive medications.
- Have any disease that is stimulated by heat, such as recurrent Herpes Simplex (within the treatment area).
- Have any active skin condition within the treatment area, such as sores, psoriasis, eczema or a rash.
- Have a history of keloid scarring, abnormal wound healing or fragile skin.
- Have had any surgical procedure, laser resurfacing or deep chemical peeling within the treatment area within the past three months or are still undergoing a process to complete healing.
- Have permanent fillers within the treatment area. If you have temporary dermal fillers within the treatment area, consult your physician before using the device.
- Have had Botox administered in the treatment area within the past three months.

- Have suspicious lesions within the treatment area.
- Currently use medications, herbal preparations, vitamins or food supplements that may cause fragile or dry skin.

POSSIBLE SIDE EFFECTS

The following side effects are rare. However, should they occur, they are usually temporary. If any of the following conditions appear and last more than 24 hours, consult your physician:

- Pain
- Damage to natural skin texture (crusting, blistering, burns, fragile skin, bruising)
- Mild headache

For other aesthetic procedures or should you have any uncertainty regarding potential side effects, consult your physician before using the device.

ABOUT TRIPOLLAR MULTI-RF TECHNOLOGY

TriPollar MULTI-RF technology delivers safe and effective skin renewal and tightening. TriPollar home use devices bring professional level aesthetic results from the professional clinic to the comfort of your home. TriPollar treatments produces real, noticeable immediate and long-lasting results. TriPollar home use devices focuses four low power beams deep into the skin to stimulate dermal activity, increasing and accelerating the natural production of collagen (the protein that is responsible for our skins' supple texture). The increase in the skin's collagen results in tighter and suppler skin.

THE TRIPOLLAR STOP MIRACLE YOUNG INTENDED USE IS AS FOLLOWS:

The TriPollar STOP MIRACLE YOUNG system is intended for non-invasive skin rejuvenation at home for delicate facial areas, specifically indicated for skin tightening in the Periorbital, Perioral and Nasolabial folds.

TRIPOLLAR STOP MIRACLE YOUNG KIT INCLUDES

1. TRIPOLLAR STOP MIRACLE YOUNG DEVICE (FIG.1)
2. TRIPOLLAR STOP MIRACLE YOUNG USB CHARGING CORD (FIG.2)
3. TRIPOLLAR PREPARATION BALM (FIG.3)
Note: After the container has been opened, the gel can be used for up to 24 months.
4. TRIPOLLAR STOP MIRACLE YOUNG USER MANUAL AND WARRANTY POLICY (FIG.4)
5. TRIPOLLAR STOP MIRACLE YOUNG QUICK START GUIDE (FIG.4)

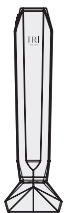
FIG.1

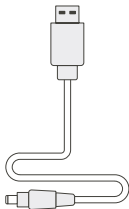
FIG.2


FIG.3

FIG.4

OVERVIEW OF THE TRIPOLLAR STOP MIRACLE YOUNG DEVICE

- | | |
|-----------------------|------------------------------------|
| 1. DEVICE BODY | 7. ON/OFF BUTTON AND INDICATOR LED |
| 2. TREATMENT HEAD | 8. GREEN INDICATOR LED |
| 3. ELECTRODE | 9. BLUE INDICATOR LED |
| 4. TEMPERATURE SENSOR | 10. PROTECTIVE COVER |
| 5. CONTROL PANEL | |
| 6. USB CONNECTION | |

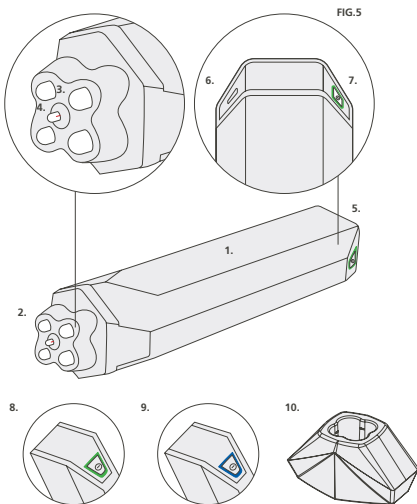THE INDICATOR LEDS ON THE TRIPOLLAR STOP MIRACLE YOUNG DEVICE WORK AS FOLLOWS:

Green indicator LED:

- Indicates that the device is ON and ready for treatment.
- Slow blinking - low battery, approximately 5 minutes remaining.
- Fast blinking - battery soon to be depleted, approximately 2 minutes remaining.

BLUE indicator LED:

- Blinking- battery charging.
- Steady- battery is fully charged.

TREATMENT PROCEDURE

COURSE OF TREATMENT

1. Treatment duration is approximately 4 minutes per area.
2. For optimal results, use your TriPollar STOP MIRACLE YOUNG device 2-3 times a week for a period of 6-8 weeks.
3. Make sure you have at least one day interval between treatments.
4. In order to maintain long term results after you have finished the full course of treatments, it is recommended that you use your TriPollar STOP MIRACLE YOUNG device 2-4 times a month or as needed.

TREATMENT PROCEDURE

1. Charge the device so that its battery is fully charged and its LED indicator is steady blue.
2. Cleanse and dry the treatment area.
3. TriPollar Preparation Balm must be used throughout treatment. Spread an even layer over the entire treatment area. In case you feel that the balm has been absorbed during the treatment, spread an additional layer of TriPollar Preparation Balm.
4. Remove the protective cover from the device. Turn TriPollar STOP MIRACLE YOUNG ON by short pressing the Power Button (FIG.6). Make sure the green indicator LED appears and the Micro Massaging effect is felt.

FIG.6

FIG.7

5. Place the device on your skin ensuring that all 4 electrodes are touching the skin.
6. Start moving the device in slow and circular motions (at a speed of about 2-4 cm/s) inside the shaded treatment area (FIG.7). You will feel pleasant warming & Micro Massaging. Make sure there is a full contact of the electrodes with the skin during treatment. Treat each area for 4 minutes.

IMPORTANT: be sure to keep the device constantly moving on your skin. Do not leave the device on one place without movement, as skin overheating may occur.

POST TREATMENT

1. Turn the device OFF by pressing the power button for two seconds.
2. Gently remove the TriPollar Preparation Balm from the skin.
3. Clean the TriPollar STOP MIRACLE YOUNG device with a soft, damp cloth or paper towel and cover it with its protective cover. Do not submerge the TriPollar STOP MIRACLE YOUNG under water!

IMPORTANT: It is required to perform a patch test 3 days prior to initial use. Apply a small amount of TriPollar Preparation Balm on the inner arm and treat the area with TriPollar STOP MIRACLE YOUNG device. Ensure that the skin does not have a negative reaction to the gel or treatment.

INDICATIONS OF SUCCESSFUL TREATMENT

As a natural result of the treatment, the skin may show some redness (erythema) and minor skin swelling (edema) in the treatment area. Both should disappear after a short time.

SAFETY MECHANISM

The device's automatic temperature-optimization mechanism protects the skin from overheating. When the skin has reached the maximal temperature, the energy source will shut down automatically.

The device automatically turns the energy source back ON, when the skin temperature drops below the maximal temperature.

TROUBLESHOOTING

In the event that your TriPollar STOP MIRACLE YOUNG device is not working, prior to contacting the service center, take the following precautions:

1. Make sure you have removed the device's protective cover before use.
2. Make sure the device's battery is charged.
3. Make sure power button was pressed and the green indicator LED appears.
4. Make sure you have applied the TriPollar Preparation Balm according to the instructions. TriPollar STOP MIRACLE YOUNG will not work properly if TriPollar Preparation Balm is not applied.

TECHNICAL SPECIFICATIONS

Technical Specifications

| | |
|--|--------------------------|
| Weight: | 60gr |
| Physical Dimensions: | H=131mm L=26mm W=26mm |
| TriPollar RF Alternating Frequencies: | 0.9 MHz / 1.1 MHz |

Operating conditions

| | |
|--------------|---------------|
| Temperature: | +5°C to +30°C |
| Humidity: | 20RH to 80RH |
| Rating: | DC 5V 2A |

Charging

When battery is fully charged (after approximately 3 hours), the LED indicator will emit a steady blue.

Transport and Storage

| | |
|--------------|--------------|
| Temperature: | 0°C to 40°C |
| Humidity: | 10RH to 90RH |

WARRANTY AND CUSTOMER SUPPORT

WARRENTY TERMS

This product has been subjected to strict quality control performed by the manufacturer. For this reason, the distributor in your country guarantees proper function of your product for 24 months from date of actual purchase. If you use the device properly, follow all instructions outlined in the user manual, you can enjoy the product to your ultimate satisfaction. This policy warrants the support that the distributor provides for the device, should it happen that failures attributable to manufacturing defects be found. In such case, you shall be entitled to make full use of the benefits of this warranty. This implies the reason and replacement of parts, free of charge (regarding both service and parts) at any Authorized Service Center. Within its validity term, you are entitled to make use of the warranty if necessary to have a fully enjoyable product, provided the above-indicated conditions are fully met. This warranty is valid for the device marketed by the distributor in your country.

The warranty is not valid when:

1. Amendments in the warranty policy, sales receipt or invoice, or the absence of any of these original documents.
2. Damage or deterioration of the device due to improper use or if the device has been modified or repaired by an unauthorized third party.
3. Damage of the product due to defects incurred during transportation, such as bumps, bangs or improper handling.
4. Damage due to humidity or any kind of liquids, foreign bodies inside the device or the use of unapproved, unadvised accessories with the device contradictory to the user manual.

5. When the device is not used in accordance with the respective user manual or with its supplied accessories.
6. Partial or total damage or loss of the device, due to natural disasters (such as earthquakes, floods, fire, lightening, etc.), force majeure or accident.
7. Damage to the device after being subjected to high dust environments, humidity or application of excessive voltage from the power inlet.
8. If the product's serial number has been amended or erased.
9. To cosmetic damage, including but not limited to scratches and dents.

CUSTOMER SUPPORT

For service, please refer to the store where you purchased the TriPollar STOP MIRACLE YOUNG or visit our website for a store locator:
www.tripollar.com | cn.pollogen.com

When reporting to local service center, please have the device serial number ready. The serial number can be found on the top panel of the device.
Customer satisfaction is our top priority. To help us provide you with the best support possible, please send us your comments and suggestions. Contact details can be found on the back cover of this user manual or via our website at: www.tripollar.com | cn.pollogen.com

APPENDIX

HAZARDOUS MATERIAL DECLARATION

TABLE 4 IS IN ACCORDANCE WITH THE PROVISIONS OF SJ/T 11364

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This table is prepared in accordance with the provisions of SJ/T 11364

- Indicates that said hazardous substance contained in all of the homogeneous materials for this part is below the limit requirement of GB/T 26572
- X Indicates that said hazardous substance contained in at least one of the homogeneous materials used for this part may be above the limit requirement of GB/T 26572

The restricted materials are used in accordance with manufacturing needs. It has no impact on device safety, and not harmful to users. The declaration is to remind end users that the device need to be disposed reasonably and environment friendly.

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© POLLOGEN LTD., ALL RIGHTS RESERVED.
TRADEMARKS OF POLLOGEN LTD.
Obelis S.A,
Obelis UK Ltd
MANUFACTURED BY POLLOGEN LTD.